

SAMARIT NURSING AIDS

For Rehabilitation and Homecare



GLIDEBOARDS
Page 2



THERAPEUTIC
BED TOPPERS
Page 4



MOVING AND HANDLING SLINGS

Page 5



TURNPLATE
TURN
CUSHION
Page 6



CAREGLIDE TRANSFER SHEETS Page 7

High quality standardsSAMARIT transfer aids have:

- Design patents
- 2 years warranty

Samarit transfer systems



SAMARIT sliding boards are Nursing aids for **moving a seated patient** with limited forces **from the middle or edge of the bed to a wheelchair, a car seat or a toilet**. Here the patient may not be able to stand but can change position when sitting independently or with the help of a healthcare professional.

While the smooth surface reduces friction resistance to a minimum, non-slip adhesive strips ensure that the GLIDEBOARD doesn't slip away. The GLIDEBOARD is made of an **unbreakable** material, is **acid-resistant** and **easy to clean** and **disinfect**. A non-slip adhesive strip enables it to be positioned safely on the bed, wheelchair, car seat or other seating options. **Five versions** are available for different purposes.

The SAMARIT GLIDEBOARD range has been developed for special use with the support of therapists.

GLIDEBOARD

SLIDING BOARD / SEATED TRANSFER



CURVED GLIDEBOARD THE GREAT THE BIGGEST SELFGLIDE TETRAGLIDE

- Unbreakable
- Sliding surfaces
- Non slide sticker
- Lightweight
- Left/right usable
- Instruction for use
- 2 years warranty

Greater mobility with GLIDEBOARD

Seated transfer for care stations, therapy centers and home care.



THE BIGGEST

Order No. 476.0102 UDI/GTIN 7640144460921



THE GREAT

Order No. 476.0080 UDI/GTIN 7640144460938



CURVED GLIDEBOARD

Order No. 476.0012 UDI/GTIN 7640144460945



SELFGLIDE

Order No. 476.0050 UDI/GTIN 7640144460952



TETRAGLIDE

Order No. 476.0863 UDI/GTIN 7640144460969



Sit-On Monofil Sliding-Towel

Order No. 476.0200 UDI/GTIN 7640144460976 This sliding board makes it easier for healthcare professionals to transfer patients from the middle of the bed to a wheelchair, to the toilet, to the car seat or a chair.

Size: 100 x 34 cm / 39 x 13 inch

SWL*: 220 kg / 485 lbs Weight: 2.9 kg / 6.4 lbs

Material: HDPE

This sliding board makes it easier for **healthcare professionals** to **transfer patients** from **edge of the bed to a wheelchair.** The product is used mainly today to transfer patients **from wheelchairs into large vehicles** (SUVs etc.).

Size: 80 x 33 cm / 32 x 13 inch

SWL*: 220 kg / 485 lbs Weight: 2.6 kg / 5.7 lbs

Material: HDPE

The original GLIDEBOARD «**Banana**» sliding board is a very tried and tested classic from Switzerland, developed with therapists by SAMARIT and patented as a design. This sliding board is often used in particular to enable **sedentary patients to transfer independently to wheelchairs, car seats or chairs**.

Size: 68 x 23 cm / 27 x 9 inch

SWL*: 200 kg / 440 lbs Weight: 1.2 kg / 2.6 lbs

Material: HDPE

This sliding board is used mainly by patients who can **transfer themselves under their own steam**.

Size: 50 x 21 cm / 20 x 8 inch

SWL*: 180 kg / 397 lbs Weight: 0.8 kg / 1.8 lbs

Material: HDPE

This sliding board has been **specially developed for paraplegic patients**. The carer can move the patient safely onto the sliding board **without lifting thanks** to its unique fish tail. It is characterised by the ease with which the patient can be sat upon, by the manner in which it can be **inserted from both sides** and is utilised for **transfer bed – chair – toilet, etc**.

Size: 63 x 23 cm / 25 x 9 inch

SWL*: 200 kg / 440 lbs Weight: 1.1 kg / 2.4 lbs

Material: HDPE

This Sit-On Monofil Sliding-Towel makes it **easy to slide onto the sliding board without clothing**. It is used by **healthcare professionals** and also by the **patients themselves**. Usually it is used with the Curved, Selfglide and Tetraglide. The monofil glide belt is suitable for **baths**, **showers** or the **toilet**.

Size: 30 x 45 cm / 12 x 18 inch

Material: Nylon fabric

Moving & Handling Slings

Handling slings are specially manufactured, upholstered cushion are for supporting a move.

Weight:



PAT-SLING

PULL-UP SLING

Order No. 450.0195 UDI/GTIN 7640144461003 In the nursing sector, the **Pat-Sling is a versitile aid** and very **easy and safe in application**. Its permit utilisation for large as well as for slim patients (up to 120 kg).

Size: 20 x 50 cm / 8 x 20 inch SWL*: 120 kg / 265 lbs

0.5 kg / 1.1 lbs

Material: Nylon fabric reinforced with double hand grips



PULL-UP 6

LIFTING AND MOVING SLING

Order No. 450.5106 UDI/GTIN 7640144461010

The Samarit Pull-Up 6 – Lifting and Moving Sling

Lifting in the bed: The patient can easily be brought onto the Pull-Up by turning the patient onto his/her side. Two people can easily carry out the moving or lifting-up of the patient without the use of much strength. The grey, slideresistant foam material serves for safe positioning of the patient.

Size: 50 x 150 cm / 20 x 59 inch

SWL*: 150 kg / 330 lbs Weight: 1 kg / 2.2 lbs

Material: Nylong fabric reinforced, hand grips (6)



SLIDE-UP WIDE

MOVING GLIDE MATTRESS

Order No. 450.5110 UDI/GTIN 7640144461027 Samarit Slide-Up wide — **Moving Glide Mattress** possesses a wider upper portion which surrounds the upper body and can be held together by a nurse with one hand. This permits the **patient to be brought safely into a seated position and moved**.

Size: 100 / 50 x 100 cm / 39 / 20 x 39 inch

SWL*: 150 kg / 330 lbs Weight: 1.6 kg / 3.5 lbs

Material: Nylon fabric with wide backpart to hold patient with one hand,

two strong handles on each side.



GLIDEBELT

Order No. 450.2100 UDI/GTIN 7640144461041

Order No. 450.2200 UDI/GTIN 7640144461058

Order No. 450.2300 UDI/GTIN 7640144461065

GLIDEBELT helps to position the slideboards under heavy patients

Small 72 - 90 cm / 28 - 35 inch

Medium 88 - 106 cm / 35 - 42 inch

Large 104 - 122 cm / 41 - 48 inch





ONE-WAY SLIDE CUSHION

Order No. 450.4060 UDI/GTIN 7640144461034 The SAMARIT Seat One-way Slide cushion only **slides in one direction**. It allows healthcare professionals to pull a patient who has slipped forward back to the correct sitting position **with little effort**.

It is ideally suited for use in **rehabilitation clinics**, care homes and care at home.

Size: 45 x 42 cm / 18 x 17 inch

Weight: 0.8 kg / 1.8 lbs

Material: Tissue with one-slide, Top PU-textile grey coating

Turnplate / Turn cushion / Transfer sheets





Nursing aids for professionals

Aid for controlled rotation, standing and seated patient transfer in clinics, nursing homes and care homes.



TURNPLATE 30
Order No. 480.0100
UDI/GTIN 7640144461089

The SAMARIT Turnplate is ideally suited to a **controlled rotation of the patient when standing up**.

Size: Ø 30 / 40 cm, Ø 12 / 16 inch, 8 mm / 0.3 inch high only

Weight: 0.5 kg / 1.1 lbs
Material: ABS / Teflon, blue



TURNPLATE 38 Order No. 480.0380 UDI/GTIN 7640144461072

The safest way to rotate a patient is when the helper puts his right foot in the middle of the Turnplate and the patient places his feet to the right and left of it.

Size: Ø 38 / 49 cm, Ø 15 / 19 inch, 8 mm / 0.3 inch high only

Weight: 0.6 kg / 1.3 lbs Material: ABS / Teflon, blue

Limited stock





BIG-TURNPLATE
Order No. 480.0400
UDI/GTIN 7640144461096

The SAMARIT BIG bariatric Turnplate is particularly suitable for a **controlled rotation of overweight patients**. The safest way to rotate a patient is when the helper puts his right foot in the middle of the turn plate and the patient places his feet to the right and left of it.

Size: Ø **41** / **54** cm, Ø 15 / 21 inch, **7** mm / 0.28 inch high only

Weight: 0.8 kg / 1.7 lbs

Material: ABS / Teflon, blue, full-scale anti-slip coating



SEAT & CHAIR TURN CUSHION

Order No. 480.1038 UDI/GTIN 7640144461102 The SAMARIT Turn Cushion is often used to **rotate a patient so they can sit on a car seat, a chair or an armchair**. It makes it easier to swing the legs round and thus makes it possible to stand up safely from the sitting position.

It is ideally suited for use in the **domestic sector**.

Size: Ø **38 cm** / Ø 15 inch **Weight: 0.12 kg** / **0.26 lbs**

Material: Antimicrobic surface, Biocompatibility DIN 10993-5, Flameretardant, blue

The SAMARIT turnplates are characterized by handiness, easy cleaning, safety and longevity.



Therapeutic Bed Toppers

The Samarit Mobicare Bed Topper has been developed to ease the **nursing of longterm bedridden patients** and is manufactured from **non-irritant** and **fire retardant materials**. **Mobilise** / **turn patient easy to avoid pressure source**.

The mobile Samarit Minicare Bed Topper cushion is small, handy and universal to use. The Minicare is used for the gentle moving and positioning in the bed. Patients who can no longer move themselves without assistance are helped by Samarit Minicare to greater mobility and independence.



MOBICARE

Mobilise / turn patient easy to avoid pressure source

Size: 90 x 195 cm / 35 x 77 inch

SWL*: 150 kg / 331 lbs

Material: with one side incontinent-proof 200 g/m2 hollowfibre

Order No. 450.0195

UDI/GTIN 7640144460990




Size: 60 x 60 cm / 24 x 24 inch breathable tissue

SWL*: 100 kg / 220 lbs Material: Microfibre yellow

Order No. 450.0260

UDI/GTIN 7640144460983

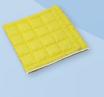

MINICARE



CareGLIDE Model «1»

(Border red)

Order No. 450.0070 UDI/GTIN 7640144461112

1

Application: seated movements to the bedside, 90° rotation, small movements, sitting up in bed.

Size: 70 x 70 cm / 28 x 28 inch

Material: 2-fold silicone coated, PE fabric, skinfriendly, top slidable, washable 74°, latex-free



CareGLIDE Model «2»

(Border yellow)


Order No. 450.0100 UDI/GTIN 7640144461126

2

Application: right/left reclining movement, upper body on the side, turning patient into side position.

Size: 100 x 80 cm / 39 x 32 inch

Material: 2-fold silicone coated, PE fabric, skinfriendly, top slidable, washable 74°, latex-free

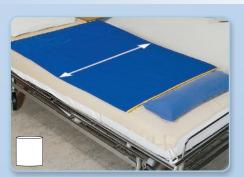

CareGLIDE Model «3»

(Border yellow)

Order No. 450.0130 UDI/GTIN 7640144461133



Application: right/left reclining movement of whole body, turning patient into side position, moving for replacing bed linen.

Size: 130 x 100 cm / 51 x 39 inch

Material: 2-fold silicone coated, PE fabric, skinfriendly, top slidable, washable 74°, latex-free



CareGLIDE Model «4»

(Border red)

Order No. 450.0180 UDI/GTIN 7640144461140



Application: right/left reclining movement, upper/lower sitting or reclining movement, transfer to bedside or wheelchair.

Size: 110 x 80 cm / 43 x 32 inch

Material: 2-fold silicone coated, PE fabric, skinfriendly, top slidable, washable 74°, latex-free



CareGLIDE Model «5»

(Border yellow)

Order No. 450.1310 UDI/GTIN 7640144461157



For heavyweight patients.

Application: Moving in bed to upper/lower reclining and sitting position, transfer bed / bed with the aid of 2 persons and pull sheet for very heavy patients, transfer / moving of spine damaged patients in bed.

Size: 130 x 100 cm / 51 x 39 inch

Material: 2-fold silicone coated, PE fabric, skinfriendly, top slidable, washable 74°, latex-free



CareGLIDE Model «6»

(Border yellow)

Order No. 450.7180 UDI/GTIN 7640144461164



For heavyweight patients.

Application: Transfer bed/bed with the aid of 2 persons and pull sheet for very heavy patients, transfer/moving of spine damaged patients in bed.

Size: 180 x 70 cm / 71 x 28 inch

Material: 2-fold silicone coated, PE fabric, skin-

Material: 2-fold silicone coated, PE fabric, skir friendly, top slidable, washable 74°, latex-free



CareGLIDE

TRANSFER SHEETS / TUBULAR SHEETS FLAT SHEETS





CareGLIDE Moving Set

Flat Sheets / Tubular Sheet Order No. 450.3160 UDI/GTIN 7640144461164

CareGLIDE Flat Sheet

Order No. 450.3162 UDI/GTIN 7640144461171 Order No. 450.3163 UDI/GTIN 7640144461188

High quality standardsSAMARIT transfer aids have:

MODEL 01 - 06 >





Moving Set consisting:


Size:

- 1 x 60 x 170 cm, 1 x 24 x 67 inch with handles
- 1 x 60 x 170 cm, 1 x 24 x 67 inch without handles
- 1 x 55 x 65 cm, 1 x 22 x 26 inch, tube

Material: 2-fold silicone coated, PE fabric, skinfriendly, top slidable, washable 74°, latex-free

Size:

- 1 x 60 x 170 cm, 1 x 24 x 67 inch with handles
- 1 x 60 x 170 cm, 1 x 24 x 67 inch without handles

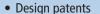

• 2 years warranty







Samarit Medical AG

Gewerbestrasse 12 CH-8132 Egg, Schweiz Tel.: +41 44 918 10 11 Fax: +41 44 918 24 39 e-mail: info@samarit.com - www.samarit.com

www.samarit.com